
GB	MASSAGER	4
RUS	МАССАЖЕР	5
CZ	MASÁŽNÍ PŘÍSTROJ	6
BG	УРЕД ЗА МАСАЖ	8
PL	MASAŻER	9
RO	DISPOZITIV PENTRU MASAJ	11
UA	МАСАЖЕР	12
SCG	MACEP	13
EST	MASSAAŽIAPARAAT	15
LV	MASĀŽAS APARĀTS	16
LT	MASSAŽIERIUS	17
H	MASSÁZSŐR	19
KZ	УҚАЛАҒЫШ	20
D	MASSAGEGERÄT	22

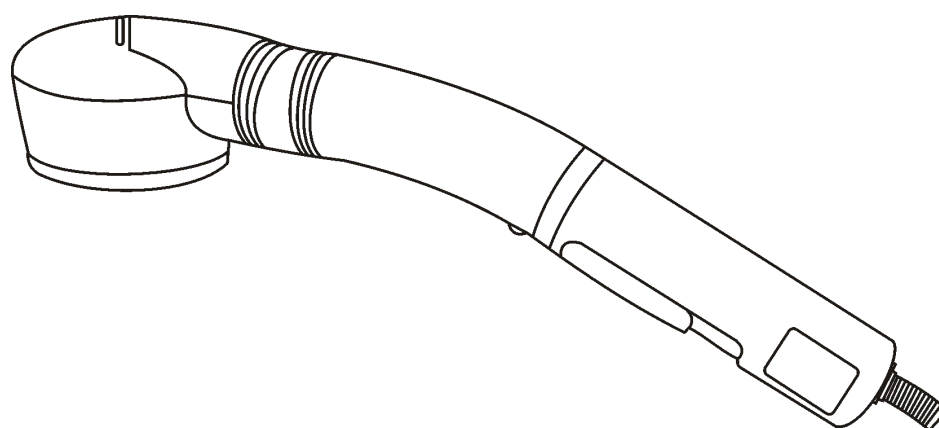GB DESCRIPTION

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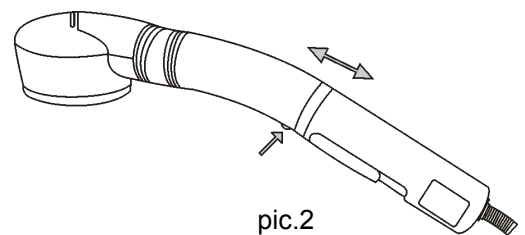
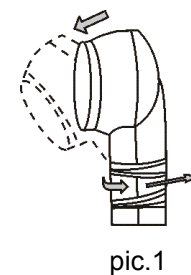
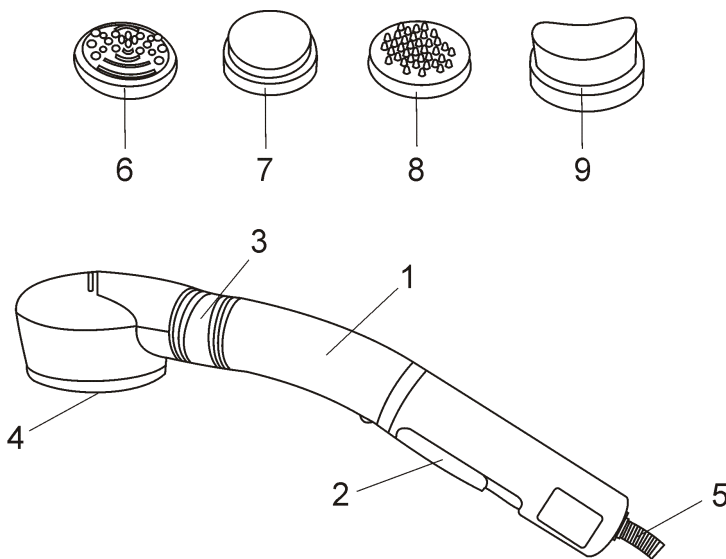
GB INSTRUCTION MANUAL

IMPORTANT SAFEGUARDS

- Please read this instruction manual carefully before use and keep in a safe place for future reference.
- Before the first connecting of the appliance check that voltage indicated on the rating label corresponds to the mains voltage in your home.
- For home use only. Do not use for industrial purposes. Do not use the appliance for any other purposes than described in this instruction manual.
- Do not use outdoors.
- Always unplug the appliance from the power supply before cleaning and when not in use.
- To prevent risk of electric shock and fire, do not immerse the appliance in water or in any other liquids. If it was happened DO NOT TOUCH the appliance, unplug it immediately and check in service center.
- Do not take this product into a bathroom or use near water.
- Close supervision is necessary when it is used near children.
- Do not leave the appliance switched on when not in use.
- Do not use other attachments than those supplied.
- Do not operate after malfunction or cord damage.
- Do not attempt to repair, adjust or replace parts in the appliance. Check and repair the malfunctioning appliance in the nearest service center only.
- Keep the cord away from sharp edges and hot surfaces.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not pull, twist, or wrap the power cord around appliance.
- Never drop or insert any object into openings on this appliance.
- Do not put operating appliance on any surface, it may cause injury or damage.
- The massager is designed for intermittent use: it can be applied without interruption not longer than **20** minutes, then the motor must be allowed to cool.
- Do not use if you suffer from thrombosis, if you pregnant or bodily injuries. Consult your doctor first.
- Immediately interrupt procedure if you feel discomfort, or if any pain or irritation results; never use for swollen or inflamed feet, or if feet have any skin eruptions. First, obligatory consult your physician.

GENERAL INSTRUCTIONS

- Use of this massager has a relaxing, soothing and stimulating effect on your muscles. If used regularly for several minutes on a daily basis, it can reduce stress, fatigue and muscular pain, and it has anti-cellulite function.
- When using skin or face remedies, this appliance helps to spread the product more evenly and thus enhance its effective skin penetration.
- This appliance is equipped with an infrared lamp. The infrared radiation penetrates the deeper layers of skin causing the pores and vessels to expand and simultaneously enhancing the metabolism process and microcirculation. Using the infrared massager brings relief from muscular pain and regenerates the skin.

ACCESSORIES

GENTLE FACIAL ATTACHMENT

- Massage face by circular motions starting from your hairline. Always avoid contact with eyes and lips. Holding taut increases the effectiveness of massage.

CONTOUR BODY ATTACHMENT

- Massage with contour body attachment improves the blood circulation and reduces accumulation of salts as well as cellulite efficiently. Use on the body directly without clothes.

MULTI-ACCUNODE MESSAGE ATTACHMENT

- This attachment provides deep and strong penetrating massage to large muscle area for relaxation and soothing muscle pain.

STIMULATING SCALP ATTACHMENT

- Scalp massage beginning at the base of the scalp slowly moving the massager up and over the crown of your head, then reverse directions for a soothing, invigorating massage.

NOTE:

- To ensure optimum effectiveness during use of any function, move the unit slowly, in either a circular or a back and forth motion.
- Always begin at a spot that is most distant from the heart. Then gradually work your way closer to the heart, using the motions described.
- Never treat the same spot longer than 3 minutes at a time.

ADJUSTMENT

ROTATION OF MESSAGE HEAD

- The massage head can be rotated downward 180° by twisting the rotation ring at the neck. Attention must be paid to twist the rotation ring in the proper direction as indicated. Otherwise it may cause damage. Do not over-twist the rotation ring (pic.1).

EXTENSION OF HANDLE

- The handle can be extended simply by pressing on the Extension control knob. And then push the handle up and down. There are 3 positions of handle extension (pic.2).

DIRECTION FOR USE

- Plug the appliance to the power supply and slide the function switch to the desired position:
- **OFF - Off.**

• **M – VIBRATION MASSAGE**

Enclosed attachments allow providing the right massage mode for each part of your body. Simply place selected attachment over the heat plate of the massager.

NOTE: Although the massager may be applied through clothing, it is most effective when acts directly on the body.

• **H – INFRA-RED HEATING**

Infrared radiation penetrates deep into muscles, enhancing blood circulation and reducing tension in muscle areas. Use the infrared heat function can also accelerate the restorative processes in the case of colds, inflammation and minor injuries.

• **HM – VIBRATION MASSAGE WITH INFRA-RED HEATING**

Combination of infrared heat and vibration massage functions is synergetic and most effective when applied directly to the skin.

CLEAN AND CARE

- Switch off the appliance and ensure that the massager is completely cool down.
- Wipe over the housing using a soft damp cloth and then - with a dry cloth.
- Do not use hard abrasives or cleaners.
- Attachments may be cleaned in soap water

STORAGE

- Be sure that the appliance is unplugged and cool down.
- Complete item CLEAN AND CARE
- Keep the appliance in a cool, dry place.

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