
ENGLISH

PROFESSIONAL BACK MASSAGER
OPERATING INSTRUCTIONS
Please read and save these instructions.
Designed by professionals - for professionals...
This device was specially designed for massaging the back and body using "chopping" movements according to the principle of karate punches...

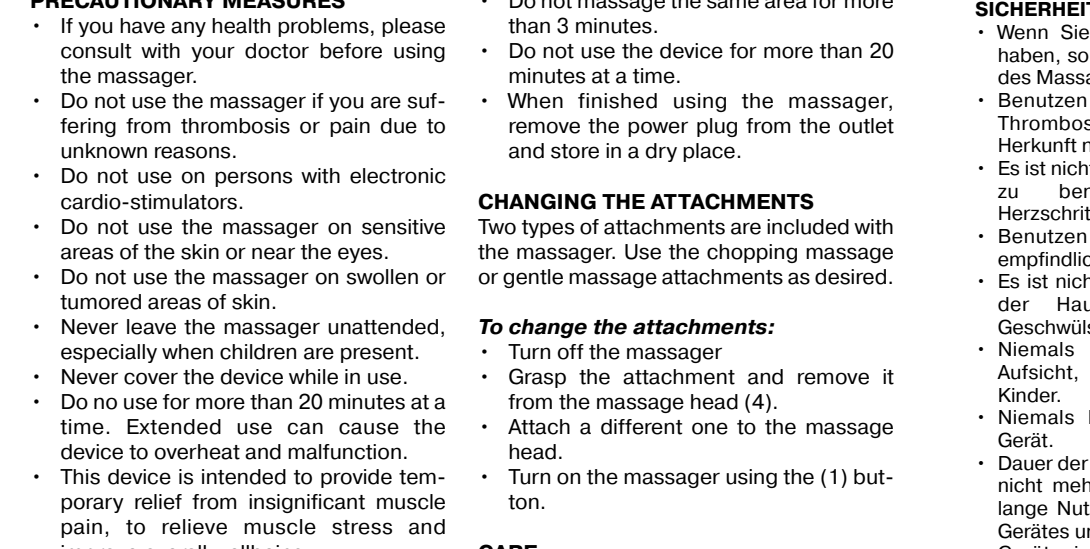PRECAUTIONARY MEASURES
If you have any health problems, please consult with your doctor before using the massager.
Do not use the massager if you are suffering from thrombosis or pain due to unknown reasons.
Do not use on persons with electronic cardio-stimulators...

CAUTIONARY MEASURES (continued)
Do not use the massager on sensitive areas of the skin or near the eyes.
Do not use the massager on swollen or tumored areas of skin.
Never leave the massager unattended, especially when children are present.

CAUTIONARY MEASURES (continued)
Do not use the massager in a shower or tub or sink.
Do not use the device in the shower or while you are bathing.
Do not use the power cord as a carrying handle for the massager.

TECHNICAL CHARACTERISTICS
Power requirements: 220-230 V 50 Hz
Power usage rating: 19 Wt
This product conforms to the EMC-Requirements as laid down by the Council Directive 89/336/EEC and to the Low Voltage Regulation (73/23/EEC)

DESCRIPTION
1. On/Off/Intensity level selector (0/1/2) Cordy
2. Body
3. Handle repositioning buttons
4. Massage heads
5. Attachments for chopping massage
6. Attachments for gentle massage
INSTRUCTIONS FOR USE
Before turning on the massager, make sure that the electrical supply voltage is compatible with the operating voltage of the massager.

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