
| | | |
|-----|------------------------------|----|
| GB | MASSAGER..... | 4 |
| RUS | МАССАЖЕР..... | 5 |
| CZ | MASÁŽNÍ PŘÍSTROJ..... | 6 |
| BG | УРЕД ЗА МАСАЖ..... | 7 |
| PL | MASAŻER..... | 8 |
| RO | DISPOZITIV PENTRU MASAJ..... | 9 |
| UA | МАСАЖЕР..... | 10 |
| SCG | MACER..... | 11 |
| EST | MASSAAŽIAPARAAT..... | 12 |
| LV | MASĀŽAS APARĀTS..... | 13 |
| LT | MASSAŽIERIUS..... | 14 |
| H | MASSÁZSŐR..... | 15 |
| KZ | УҚАЛАҒЫШ..... | 16 |
| D | MASSAGEGERÄT..... | 18 |
| CR | APARAT ZA MASAŽU..... | 19 |

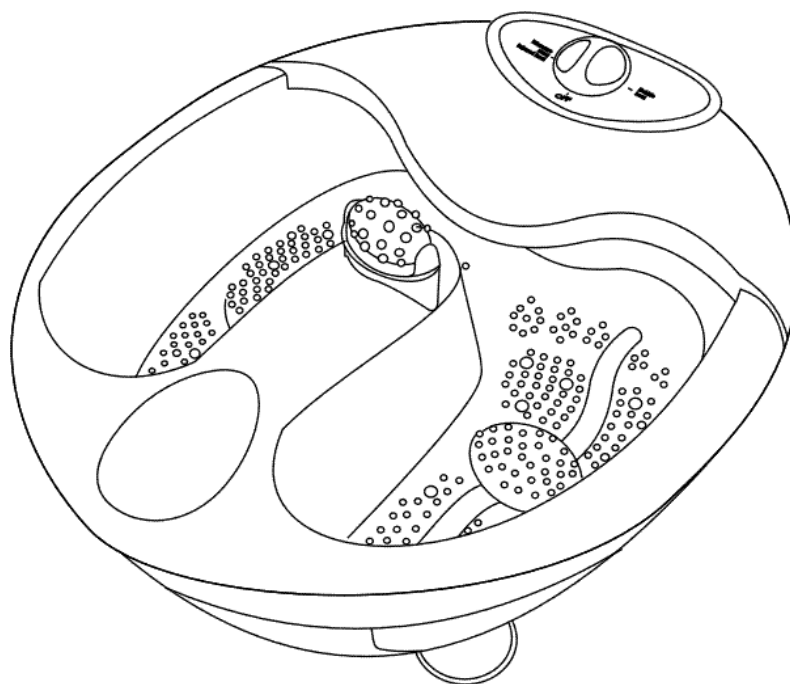GB DESCRIPTION

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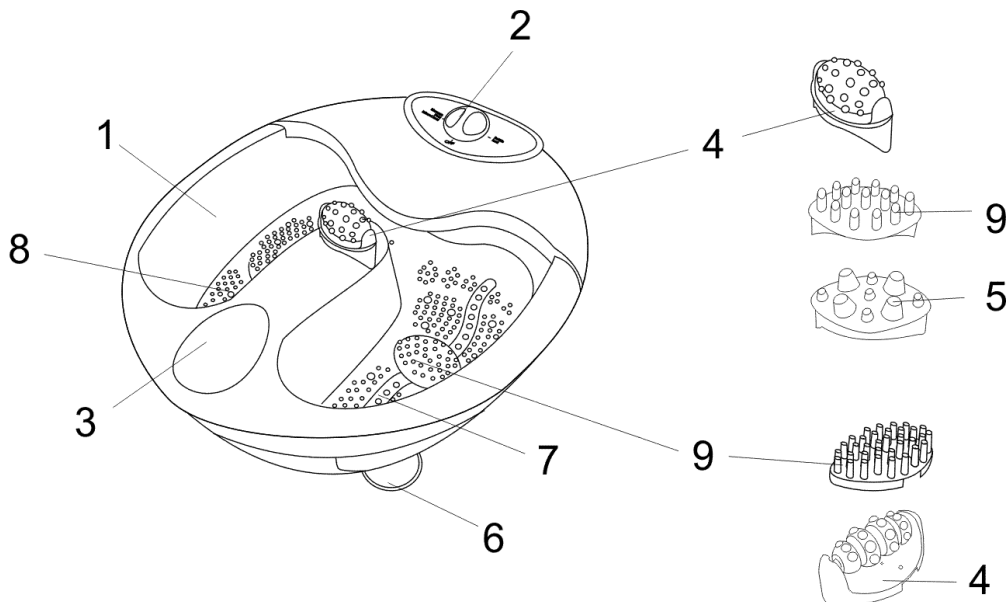
GB INSTRUCTION MANUAL

IMPORTANT SAFEGUARDS

- Please read this instruction manual carefully before use and keep in a safe place for future reference.
- Before the first connecting of the appliance check that voltage indicated on the rating label corresponds to the mains voltage in your home.
- For home use only. Do not use for industrial purposes. Do not use the appliance for any other purposes than described in this instruction manual.
- Do not use outdoors.
- Always unplug the appliance from the power supply before cleaning and when not in use.
- To prevent risk of electric shock and fire, do not immerse the appliance in water or in any other liquids. If it was happened DO NOT TOUCH the appliance, unplug it immediately and check in service center.
- Do not take this product into a bathroom or use near water.
- Close supervision is necessary when it is used near children.
- Do not leave the appliance switched on when not in use.
- Do not use other attachments than those supplied.
- Do not operate after malfunction or cord damage.
- Do not attempt to repair, adjust or replace parts in the appliance. Check and repair the malfunctioning appliance in the nearest service center only.
- Keep the cord away from sharp edges and hot surfaces.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not pull, twist, or wrap the power cord around appliance.
- The massager is designed for intermittent use: it can be applied without interruption not longer than **20** minutes, then the motor must be allowed to cool.
- Do not use the massager if you are pregnant, or suffer from thrombosis, open wounds, broken skin and varicose veins. Consult your doctor first.
- Do not attempt to plug in or unplug the unit while feet are in water.
- Do not stand up or place any objects in the unit, use it only when you are sitting.
- Immediately interrupt procedure if you feel discomfort, or if any pain or irritation results; never use for swollen or inflamed feet, or if feet have any skin eruptions. First, obligatory consult your physician.

GENERAL INSTRUCTIONS

- Foot massager bath at regular intervals can improve blood circulation in the feet and indirectly in the entire leg area as well.
- Temperature-controlled full-length heaters under each of footrests maintain the desired temperature.
- There are **40** air jets provide a relaxing warm bubble bath for your feet.
- This appliance is equipped with an infrared lamp. The infrared radiation penetrates into the deep layers of skin causing the pores and vessels to expand and simultaneously enhancing the metabolism process and microcirculation. Using the infrared massager brings relief from muscular pain and regenerates the skin.

DIRECTION FOR USE

- Place the appliance on a level surface and the rubber feet will prevent it from slipping.
- Be sure the appliance is switched off (**OFF**).
- Pour warm water (in desired temperature) into the unit. Do not pour water over the "max" water-level mark. This unit is designed to maintain the water temperature, not to raise it.

DO NOT PLUG THE APPLIANCE IN WHILE YOUR FEET ARE INSIDE THE UNIT, WHENEVER YOU ARE IN DRY OR WET OPERATION.

- For warming-up muscles you can use the infra-red lamp located in the centre of massager, set the mode switch to the according position.
- Place feet on the foot rests and select your desired setting:

WARMTH AND MASSAGE + INFRARED HEATING

- Set control switch to (**Massage / Heat / Infrared**) mark if you want to massage your feet in warm water. A foot vibration massage without water is also possible with this setting.
- Select (**Massage / Bubbles / Heat / Infrared**) setting if you want to have the unit filled with warm water for a foot massage and bubble bath. The built-in heater delays temperature decrease of the warm water.

WARMTH

- Warm water footbath with bubbling is ideal for relaxation after previous massage. Set switch to (**Bubbles / Heat**).

REMEMBER NOT TO STAND ON THE UNIT.

- When you are through with the massager, just switch it off. Then remove your feet and unplug the unit.

CLEAN AND CARE

- Switch off the appliance and ensure that the massager is completely cool.
- Do not use hard abrasives or cleaners.
- Attachments may be cleaned in soapy water.
- After using, simply rinse out with water. Never immerse the entire unit into water. If you wish, you can wipe off all surfaces with a cloth damped with water and a mild detergent, then rinse and dry.

STORAGE

- Be sure that the appliance is unplugged and cool down.

- Complete item CLEAN AND CARE.
- Keep the appliance in a cool, dry place.

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