
GB	MASSAGER.....	4
RUS	МАССАЖЕР.....	5
CZ	MASÁŽNÍ PŘÍSTROJ.....	7
BG	УРЕД ЗА МАСАЖ.....	8
PL	MASAŻER.....	10
RO	DISPOZITIV PENTRU MASAJ.....	11
UA	МАСАЖЕР.....	13
SCG	MACEP.....	14
EST	MASSAAŽIAPARAAT.....	16
LV	MASĀŽAS APARĀTS.....	17
LT	MASSAŽIERIUS.....	18
H	MASSÁZSŐR.....	20
KZ	УҚАЛАҒЫШ.....	21
D	MASSAGEGERÄT.....	23

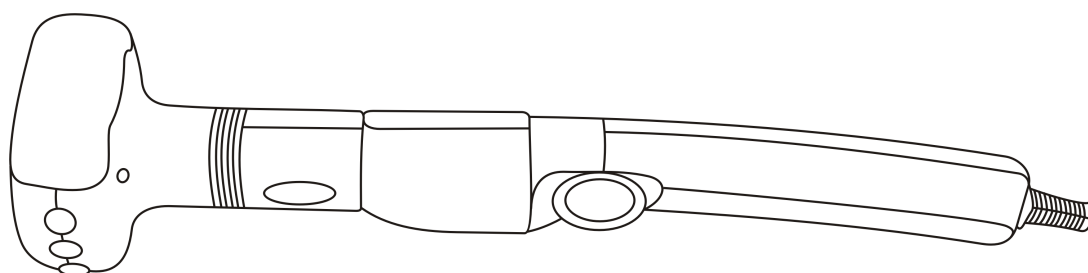GB DESCRIPTION

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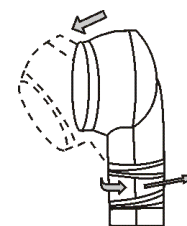
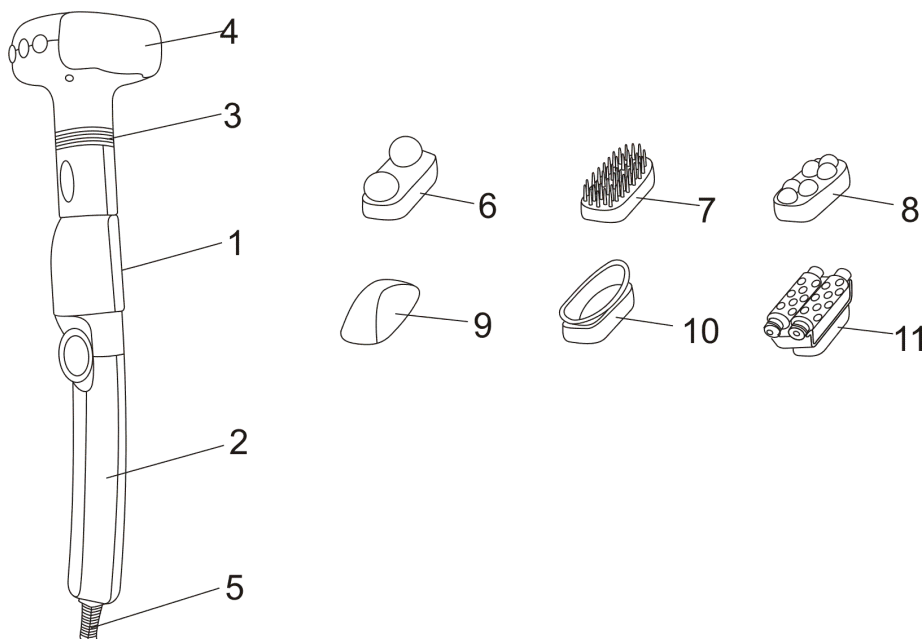
ADJUSTMENT

ROTATION OF MESSAGE HEAD

- The massage head can be rotated downward 270° by twisting the rotation ring at the neck. Attention must be paid to twist the rotation ring in the proper direction as indicated. Otherwise it may cause damage. Do not over-twist the rotation ring (pic.1).

TWISTED HANDLE

- You can turn the handle on 120° for simplification of access to any sites of a body

DIRECTION FOR USE

- Plug the appliance to the power supply and slide the function switch to the desired position:
- **0 - Off.**

- **1 – VIBRATION MASSAGE**

Enclosed attachments allow providing the right massage mode for each part of your body. Simply place selected attachment over the heat plate of the massager.

NOTE: Although the massager may be applied through clothing, it is most effective when acts directly on the body.

- **2 – INFRA-RED HEATING**

Infrared radiation penetrates deep into muscles, enhancing blood circulation and reducing tension in muscle areas. Use the infrared heat function can also accelerate the restorative processes in the case of colds, inflammation and minor injuries.

- **3 – VIBRATION MASSAGE WITH INFRA-RED HEATING**

Combination of infrared heat and vibration massage functions is synergetic and most effective when applied directly to the skin.

CLEAN AND CARE

- Switch off the appliance and ensure that the massager is completely cool down.
- Wipe over the housing using a soft damp cloth and then - with a dry cloth.
- Do not use hard abrasives or cleaners.
- Attachments may be cleaned in soap water

STORAGE

- Be sure that the appliance is unplugged and cool down.
- Complete item CLEAN AND CARE
- Keep the appliance in a cool, dry place.

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