
Open your package

- Prestigio Smart Body Mass Scale
- Four carpet feet
- Warranty card
- Four AAA-size carbon batteries (1.5V each)
- Quick start guide
- Legal and safety notice

Device overview

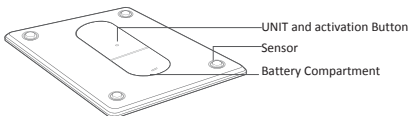LCD display

kg Kilogram  Data transmitting

lb Pound  Data pending to transmit to wireless wellness system

st Stone  Low Battery

Getting started

1. Insert the Batteries

1. Open the battery door in the back of the scale.
2. Insert the batteries (4 x 1.5V AAA) into the battery compartment according to the polarity indications marked inside the compartment.
3. Close the battery door and wait until the digits "0.0kg" is shown on the LCD.

2. Prestigio Family Health Compass application

You may store, monitor and share measurement data on your mobile device with Prestigio Family Health Compass application. To learn more about this application and download the latest version, scan the QR code at right with your mobile device to visit:

http://www.prestigio.com/hc_app_redirect.html

Note: iOS application is compatible with Apple devices equipped with Bluetooth SMART technology (iPhone 4S and newer, iPad generation 3 and newer, iPad Mini, iPod generation 5 and newer).

3. Pair-up the Scale with Mobile Device

1. Turn on Bluetooth and App. Make sure both are ON when pair-up is proceeding.
2. Press and hold "UNIT" button in the back of the scale to start pair-up.
 - Symbol  and symbol  will be shown on the LCD alternatively, indicating pair-up is proceeding.
 - If SUCCEED, symbol  will be shown on the LCD.
 - If FAIL, symbol "E1" will be shown on the LCD.

Note: The scale will return to zero, whether the pair-up is successful or failed, it will turn off automatically after steadying 60 seconds.

Measurement

4. Taking Measurements on Carpet

For using Prestigio Smart Body Mass Scale, please add the

carpet feet as demonstrated in the following page:

EN

1. Remove the non-slip pads from the scales feet.
 2. Fit the carpet feet firmly onto the scale's feet.
- Note: For the first removal of non-slip pads, you may need to use the assistance of a knife to make a small cut between the edge of the scale feet and the non-slip pad.

5. Select Measurement Unit

With batteries correctly installed, press "UNIT" button in the back of the scale to select measurement unit. The default measurement unit is "kilogram". You may press "UNIT" button to choose among kilogram, stone and pound.

6. Daily Measurement

1. With original SENSE ON patent technology, the scale will automatically switch on as you step on the platform barefooted.
2. Stand still and keep full contact with the scale until the LCD displays the weight and stops blinking the unit.
3. When your scale is successfully paired with your mobile device and the Bluetooth is ON, it will process data transmission automatically.

Data management

7. Data Transmission

With Smart Body Mass Scale successfully pair-up with your mobile device, the measurement data will be automatically transmitted to your mobile device via Bluetooth.

- The symbol  will disappear after successful data transmission, and you may check your personal health data stored in your mobile device.
- If the data transmission fails, the symbol  will remain. The pending measurement data will be temporarily kept in

the scale and transmitted to your mobile device when next measurement is completed.

Technical Specifications

Item	Specification
Dimension	Scale: 321 x 321 x 23.5 mm, Panel: 321 x 321 x 6 mm
Net Weight	Approximately 1.98 kg (Excluding the dry cells)
Display	Blue LCD with White Backlight V.A.: 74 x 53 mm
Measurement Unit	Kilogram / Stone / Pound
Measurement Range	150 kg / 330 lb
Division	0.1 kg / 0.2 lb
Accuracy	0-50 kg: ± 0.3 kg; 50-100 kg: ± 0.4 kg; 100-150 kg: ± 0.5 kg;
Working Environment	Temperature: 0 °C to 40 °C Humidity: ≤ 90 % RH
Storage Environment	Temperature: -20 °C to 60 °C Humidity: 10 % RH to 93 % RH
Power Source	6V (Four AAA-size carbon batteries)
Auto-ON	SENSE ON technology
Auto-OFF	About 7-13 seconds while showing 0.0 About 10-20 seconds after the weight data is locked
Protection Against Ingress of Water	IPX0
Software Version	1.0

Disclaimer

CE As the Prestigio products are constantly updated and improved, your device's software and hardware may have a slightly different appearance or modified functionality than presented in this Quick Start Guide.

EN



Prestigio hereby states that this device complies to basic requirements and another relevant regulations given by the EU directive 93/42/EEC. Full version of the declaration of conformity is available on:

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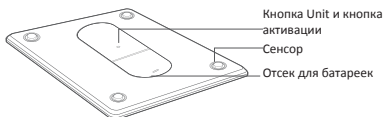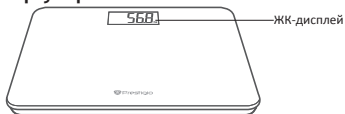2. Приложение Prestigio Family Health Compass application

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