
Open your package

- Prestigio Smart Body Fat Scale
- Four AAA-size carbon batteries (1.5V each)
- Four carpet feet
- Quick start guide
- Warranty card
- Legal and safety notice

Device overview

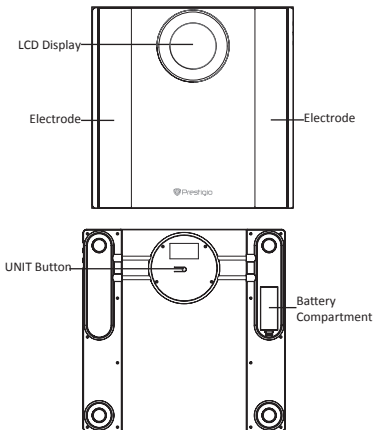LCD display

kg Kilogram

Successful Bluetooth Connection

stlb Stone

Data transmitting / pending to transmit to wireless wellness system

lb Pound

Low Battery

Getting started

1. General Instructions

Prestigio Smart Body Fat Scale lets a small amount of weak (less than 1mA) current flow through the human body, so as to detect the bio-impedance and estimate body fat.

However, please be aware that anyone with an wearable or implantable medical electronic instrument, such as a pacemaker, must avoid using this device.

2. Insert the Batteries

- 1) Open the battery door in the back of the scale.
- 2) Insert the batteries (4 x 1.5V AAA) into the battery compartment according to the polarity indications marked inside the compartment.
- 3) Close the battery door and wait until the digits "0.0kg" is shown on the LCD.

3. Prestigio Family Health Compass application

You may store, monitor and share measurement data on your smartphone or tablet with Prestigio Family Health Compass application. To learn more about this application and download the latest version, scan the QR code at right with your mobile device to visit:

http://www.prestigio.com/hc_app_redirect.html

Note: iOS application is compatible with Apple devices equipped with Bluetooth® SMART technology (iPhone 4S and newer, iPad generation 3 and newer, iPad Mini, iPod generation 5 and newer).

Android Application is compatible with devices equipped with Bluetooth® SMART technology and running Android version 4.3 or higher.

Please see <http://www.prestigio.com/healthcare/compatibility> for details.

4. Pair-up the Scale with Mobile Device

- 1) Turn on Bluetooth and App. Make sure both are ON when pair-up is proceeding.

2) Press and hold "UNIT" button in the back of the scale to start pair-up.

- Symbol  and symbol  will be shown on the LCD alternatively, indicating pair-up is proceeding.
- If SUCCEED, symbol [] will be shown on the LCD.
- If FAIL, symbol "E1" will be shown on the LCD.

Measurement

5. Taking Measurements on a Carpet

For using Prestigio Smart Body Fat Scale, please add the carpet feet as demonstrated in the following page:

- 1) Remove the non-slip pads from the scales feet.
- 2) Fit the carpet feet firmly onto the scale's feet.

Note: For the first removal of non-slip pads, you may need to use the assistance of a knife to make a small cut between the edge of the scale feet and the non-slip pad.

6. Select Measurement Unit

With batteries correctly installed, press "UNIT" button in the back of the scale to select measurement unit. The default measurement unit is "kilogram". You may press "UNIT" button to choose among kilogram, stone and pound.

7. Daily Measurement

- 1) With original SENSE ON patent technology, the scale will automatically switch on as you step on the platform barefooted.
- 2) Stand still and keep full contact with the electrodes until the LCD stops displaying moving "0".
- 3) When your scale is successfully paired with your mobile device and the Bluetooth is ON, it will process data transmission automatically.

Data management

8. Data Transmission

With Smart Body Fat Scale successfully pair-up with your mobile device, the measurement data will be automatically transmitted to your mobile device via Bluetooth.

- The symbol  will disappear after successful data transmission, and you may check your personal health data stored in your mobile device.
- If the data transmission fails, the symbol  will remain. The pending measurement data will be temporarily kept in the scale and transmitted to your mobile device when next measurement is completed.

Technical Specifications

Item	Specification
Dimension	Scale: 310 x 310 x 24.2 mm, Panel: 310 x 290.5 x 6 mm
Net Weight	Approximately 1.8 kg (Excluding the dry cells)
Display	Digital LCD with White Backlight V.A.: Round, 70 mm in diameter
Measurement Unit	Kilogram / Stone / Pound
Measurement Range	5 kg to 180 kg / 0 st: 11 lb to 28 st: 5 lb / 11 lb to 397 lb
Division	0.1 kg / 0.2 lb
Accuracy	0-50 kg: ± 0.3 kg; 50-100 kg: ± 0.4 kg; 100-150 kg: ± 0.5 kg; 150-180 kg: ± 0.7 kg
Working Environment	Temperature: 0 °C to 40 °C Humidity: ≤ 90 % RH
Storage Environment	Temperature: -20 °C to 60 °C Humidity: 10 % RH to 93 % RH
Power Source	6V (Four AAA-size carbon batteries)
Auto-ON	SENSE ON technology
Auto-OFF	About 10 seconds while showing 0.0 About 15 seconds after the weight data is locked
Protection Against Ingress of Water	IPX0

Disclaimer

Please be aware that this device is designed for adults' self-measuring and self-monitoring body fat level. Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening. This device should not be used by anyone who is acutely or chronically ill, suffering from a disease or taking medications that affect your water levels.

 As Prestigio products are constantly updated and improved, your device's software and hardware may have a slightly different appearance or modified functionalities than presented in this quick start guide.

Prestigio hereby states that this device complies to basic requirements and other relevant regulations given by the EU directive 93/42/EEC. Full version of the declaration of conformity is available on:

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