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
www.vitek.ru

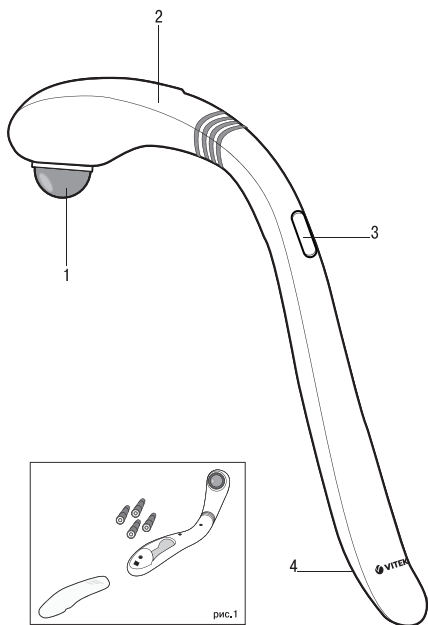HANDHELD MASSAGER VT-1393 W

The massager is intended for the collar zone massage.

DESCRIPTION

1. Massage head
2. Body
3. Power switch/operation mode switch «ON/OFF»
4. Battery compartment lid

ATTENTION!

Before using the massager consult your doctor, if you have:

- cardiovascular diseases;
- an installed pacemaker;
- oncological diseases;
- diabetes;
- spine implants;
- active chronic diseases;
- damaged skin or a skin disease in the massage area;

and if you:

- are pregnant;
- have recently undergone a surgery;
- are under treatment;
- have physical activity limits.

SAFETY MEASURES

Before using the massager, read this instruction manual carefully and keep it for future reference.

Use the unit for its intended purpose only, as specified in this manual. Mishandling the unit can lead to its breakage and cause harm to the user or damage to his/her property.

- Continuous operation time while using the massager should not exceed 10-15 minutes; make at least a 30-minute break between massage sessions.
- Do not use the unit near heat sources or open flame.
- Do not immerse the massager into water or any other liquids.
- Do not use the massager in rooms with high temperature and high humidity.
- Do not use the massager outdoors.
- Do not use the massager when you are sick or drowsy.
- The massage should give only pleasant sensations. If you feel pain or sickness during the massage, stop using the massager immediately and consult your doctor.
- Do not leave the operating unit unattended.
- Always switch the unit off before cleaning or when you are not using it.

ENGLISH

- For children safety reasons do not leave polyethylene bags used as a packaging unattended.

Attention! Do not allow children to play with polyethylene bags or packaging film. **Danger of suffocation!**

- Do not allow children to use the unit as a toy.
- During the unit's operation and breaks between operation cycles, place the unit out of reach of children.
- The unit is not intended for usage by physically or mentally disabled persons (including children) or by persons lacking experience or knowledge if they are not under supervision of a person who is responsible for their safety or if they are not instructed by this person on the usage of the unit.
- Never use the unit if it works improperly or after it was dropped.
- Do not attempt to repair the unit. Do not disassemble the unit by yourself. If any malfunction is detected or after it was dropped, apply to any authorized service center at the contact addresses given in the warranty certificate and on the website www.vitek.ru.
- To avoid damages, transport the unit in the original package only.
- Keep the unit in a dry cool place away from children and disabled persons.

THE UNIT IS INTENDED FOR HOUSEHOLD USAGE ONLY

BEFORE THE FIRST USE

After unit transportation or storage at low temperature it is necessary to keep it for at least 3 hours at room temperature before switching on.

- Unpack the unit and remove any stickers that can prevent its operation.
- Check the unit for damages, do not use it in case of damages.
- Wipe the massager body (2) with a soft dry cloth.

Battery installation (pic. 1)

- Open the battery compartment lid (4) and insert 4 AA batteries (not supplied), strictly following the polarity.
- Close the battery compartment lid (4).

Note:

- do not use old and new batteries together or different types of batteries;
- when the massager operation is unstable, check the batteries and replace them if necessary;
- remove the batteries if you are not planning to use the massager for a long time.

ENGLISH

The massager is intended for using on the collar zone.

The collar zone is easy to damage by pinching a nerve, straining a ligament etc. Besides the pharmaceutical treatment you can perform massage sessions with the massager.

MASSAGING

The massager is intended for the collar zone massage.

Besides the pharmaceutical treatment you can perform massage sessions with the massager.

Attention! *Avoid direct impact on the vertebrae and perform all massage techniques on both sides of the spine during spine massage.*

- Switch the massager on with the switch (3) and select the required massage intensity.
- Massage the selected body area no longer than 3 minutes, then move the massager to another body area.
- After a 10-15 minute massage session switch the unit off by setting the switch (3) to the position «OFF».

Note: *Make at least a 30-minute break after every massage session.*

RECOMMENDATIONS

- Use the massager in at least an hour after meals.
- Do not use the unit for head massage.
- Do not use the massager at ambient temperature above 35°C and in rooms with high humidity.
- If you feel unpleasant symptoms during using the massager, stop using the unit immediately and see your doctor.
- For massaging put the massager to the wanted body area, do not press on the massage head (1) with excessive effort.
- It is recommended to start massaging from the neck, then moving to the back.
- Recommended massage time for one body area is from 3 to 5 minutes.
- General massage session should be no longer than 10-15 minutes with a subsequent 30-minute break.

CLEANING AND CARE

- Switch the unit off by setting the power switch (3) to the position «OFF».
- Wipe the massager body (2) with a soft dry cloth; if it is very dirty, clean the unit body (2) with a slightly damp cloth and then wipe it dry.
- Do not use abrasive substances and solvents to clean the unit body (2).

ENGLISH

- Do not immerse the massager into water or other liquids.
- Remove the batteries from the battery compartment if you are not planning to use the massager for a long time.
- Keep the unit away from children in a dry cool place.

DELIVERY SET

Massager – 1 pc.

Instruction manual – 1 pc.

TECHNICAL SPECIFICATIONS

DC supply voltage: 6 V – 4 AA batteries (not supplied).

Continuous operation time: not longer than 15 minutes.

The manufacturer preserves the right to change the specifications of the unit without a preliminary notification.

Unit operating life is 3 years.

RECYCLING

For environment protection do not throw out the unit and the batteries

with usual household waste after its service life expiration; apply to the specialized center for further recycling.

The waste generated during the disposal of the unit is subject to mandatory collection and consequent disposal in the prescribed manner.

For further information about recycling of this product apply to a local municipal administration, a disposal service or to the shop where you purchased this product.

Guarantee

Details regarding guarantee conditions can be obtained from the dealer from whom the appliance was purchased. The bill of sale or receipt must be produced when making any claim under the terms of this guarantee.

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