**HOMEDICS®**

# Programmable Massager with Heat

**Instruction Manual**

---

BKP-300-2EU  
BKP-300-2GB

**2** year  
guarantee

READ ALL INSTRUCTIONS BEFORE USE.  
SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

## IMPORTANT SAFEGUARDS

---

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Keep cord away from heated surfaces.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- DO NOT use outdoors.
- DO NOT crush. Avoid sharp folds.
- This product requires a 230V AC power supply.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. For service, send to a HoMedics Service Centre. All servicing of this appliance must be performed by authorised HoMedics service personnel only.

## SAFETY CAUTIONS:


---

PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.

- **If you have any concerns regarding your health, consult a doctor before using this appliance.**
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.

- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- If you feel any discomfort whilst using this appliance, discontinue use and consult your doctor.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time.

## PRODUCT FEATURES:
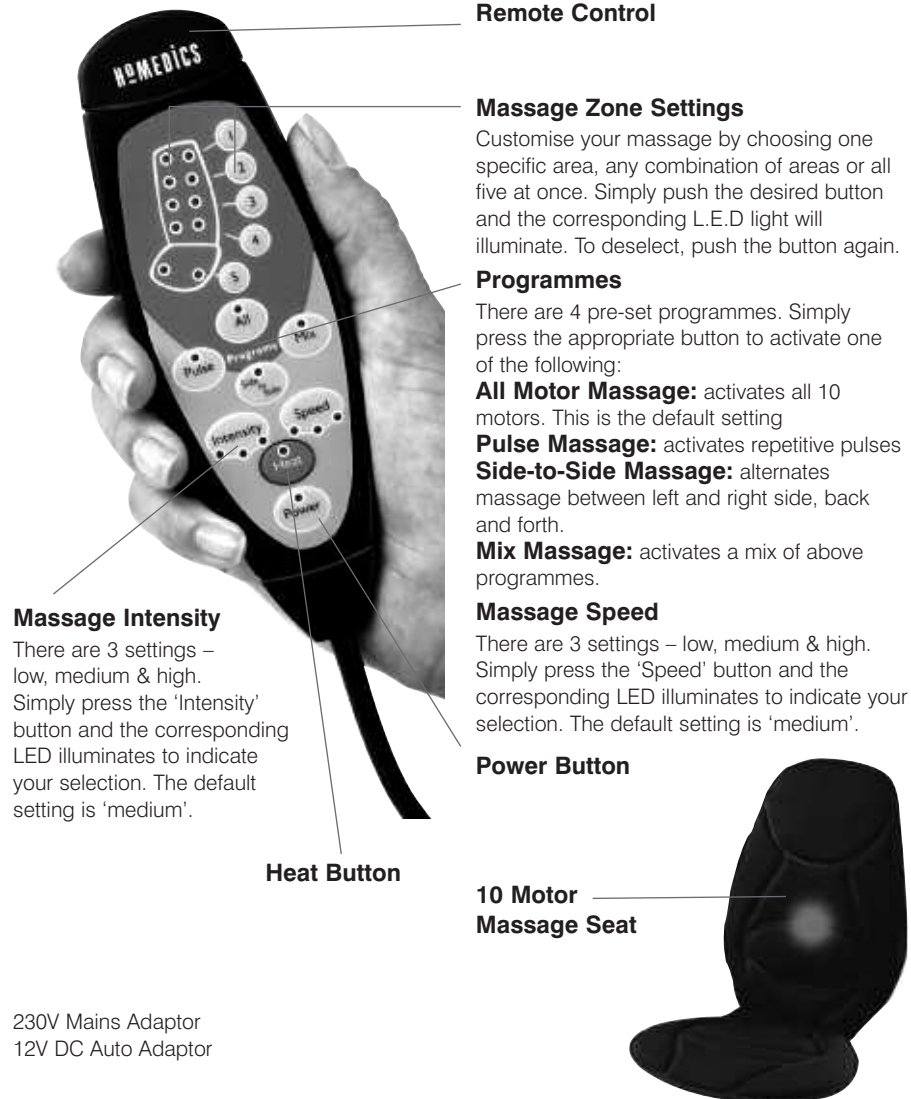
### Remote Control

### Massage Zone Settings

Customise your massage by choosing one specific area, any combination of areas or all five at once. Simply push the desired button and the corresponding L.E.D light will illuminate. To deselect, push the button again.

### Programmes

There are 4 pre-set programmes. Simply press the appropriate button to activate one of the following:

**All Motor Massage:** activates all 10 motors. This is the default setting

**Pulse Massage:** activates repetitive pulses

**Side-to-Side Massage:** alternates massage between left and right side, back and forth.

**Mix Massage:** activates a mix of above programmes.

### Massage Speed

There are 3 settings – low, medium & high. Simply press the 'Speed' button and the corresponding LED illuminates to indicate your selection. The default setting is 'medium'.

### Power Button

### 10 Motor Massage Seat

### Massage Intensity

There are 3 settings – low, medium & high. Simply press the 'Intensity' button and the corresponding LED illuminates to indicate your selection. The default setting is 'medium'.

### Heat Button

230V Mains Adaptor  
12V DC Auto Adaptor

## INSTRUCTIONS FOR USE:

### For use in the office or home

NOTE: complete steps 1 to 3 before switching the appliance on at the mains.

1. Attach the Massage cushion to almost any chair using the integrated strap located at the back of the seat. Ensure it is held firmly in place by adjusting the strap as necessary.
2. Connect the power supply lead from the adaptor with the corresponding lead in the side of the seat.
3. Plug the adaptor into a 230V mains outlet and switch on.
4. Once seated, use the remote control to operate the appliance. Press the 'Power' button and select the desired massage zone and heat button (if desired) to start the massage. Refer to the illustration of the remote control for further explanation on the different programmes available.
5. Press the 'Power' button to stop the massage.

### For use in the car

1. Attach the Massage seat using the integrated strap located at the back of the seat. Ensure it is held firmly in place by adjusting the strap as necessary.
2. Connect the power supply lead from the 12V auto adaptor with the corresponding lead in the side of the seat.
3. Plug the adaptor into a 12V power supply in the car.
4. Refer to steps 4 & 5 in the section above for further instructions for use.

NOTE:

- The 'Heat' function is pressure activated and will not function unless pressure is applied to the base of the seat. For safety purposes, heat will only function when massage is activated.
- Only gentle force should be exerted against the seat in order to eliminate risk of injury.
- There is a 15 minute auto shut-off on this appliance for your safety. This feature should not be considered as a substitute for 'OFF'. Always remember to switch the appliance off when not in use.
- Take care to position the product where it will not damage wooden or upholstered surfaces.

## MAINTENANCE:

### Cleaning

Unplug the appliance and allow to cool before cleaning. Clean only with a soft, slightly damp sponge.

- Never allow water or any other liquids to come into contact with the appliance.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

### Storage

Place the appliance in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

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