# DUO *sport*

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*User Manual - Manual de usuario - Benutzerhandbuch - Gebruiks Handleiding  
Manual do Utilizador - Manuale d'istruzioni - Instrukcja obsługi - Használati  
utasítás - Användarmanual - Ohjekirja - Brugermanual - Brukermanual -  
Руководство пользователя - دليل المستخدم - מדרוך למשתמש - دليل المستخدم*

Français .....	5
English .....	12
Español .....	19
Deutsch .....	26
Nederlands .....	33
Português .....	40
Italiano .....	47
Polski .....	54

Magyar .....	61
Svenska .....	68
Dansk .....	75
Suomi .....	82
Norsk .....	89
Русский .....	96
ةيبرعلا .....	103
תירובע .....	110
TENS instructions .....	117

# Merci !

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# Thank you !

We thank you for purchasing a Duo Sport device.

Duo Sport is a medical device, to make best use of it, you will find all relevant informations included in this user manual.

For any additional information or to simply give us your feedback, you can contact us at [contact@bluetens.com](mailto:contact@bluetens.com).

Getting started with Duo Sport is easy, here is some advice to let you benefit from it as quickly as possible and safely.

1. Please read the safety instructions as well as the contraindications at page 17. ⚠

2. Your Smartphone should support the minimum system requirements :

**iOS** : iPhone 5s or above, iOS 7 or above

**Android** : Bluetooth 4.0 compatible and Android 6 or above

3. You will need to install the Bluetens App on your Smartphone.

This application is available under « Bluetens Get Better »

on Appstore & Google play.

4. Your Bluetooth must be ON.



It is not necessary to pair your Duo Sport, your Smartphone will connect automatically as you will launch the program. This means that you will not need to fill in a PIN code or change your Bluetooth settings, Duo Sport is small, strong and clever too!

**Are all steps validated? So now we are ready to go through more serious things ;)**

## THE DUO SPORT

- 1 ON/OFF Button
- 2 Intensity adjustment knob
- 3 Program intensity indicators
- 4 Charge indicator
- 5 Power on indicator
- 6 Sport clip plug
- 8 Micro USB plug
- 9 Reset button
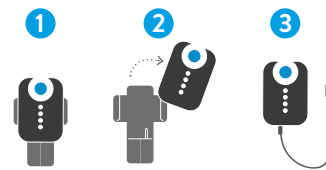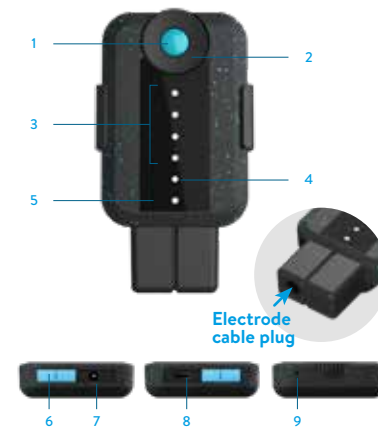
### To charge your Duo Sport :

1. Unclip the Duo sport from the main electrode
2. Remove the device from the clip
3. Push the slider to the right and connect the usb cable to charge the device


## THE APPLICATION

- AUTOMATIC PAIRING
- FAVOURITES
- DOWNLOAD

- SPLIT
- SKIP
- LINK/UNLINK


## FIRST USE
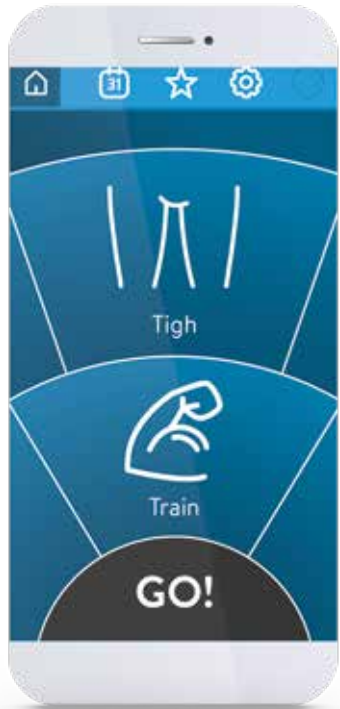
## LAUNCH THE BLUETENS APPLICATION

- > Select your device Duo Sport
  - > Creation of an account
  - > A tutorial will be available during your first use
  - > Program selection
- >Rotate the wheel to select the desired body part, then do the same for the desired action.
- + Cure**  
For physiotherapy care aiming to amplify rehabilitation or allow a rapid pain relief effect.
  - + Relax**  
For a relaxing massage or to ease muscle tension. It is also very effective to improve the quality and speed of recovery.
  - + Tone**  
For a muscular work-out session.
- >Then press «GO».

## PROGRAM DESCRIPTION

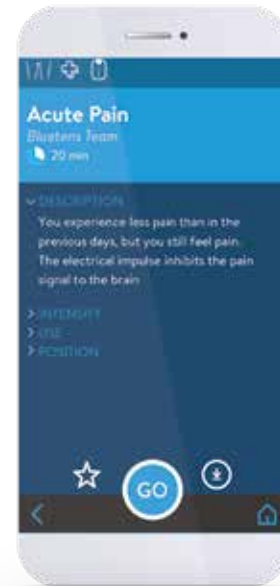
## ELECTRODE PLACEMENT

## TURN YOUR DUO SPORT ON
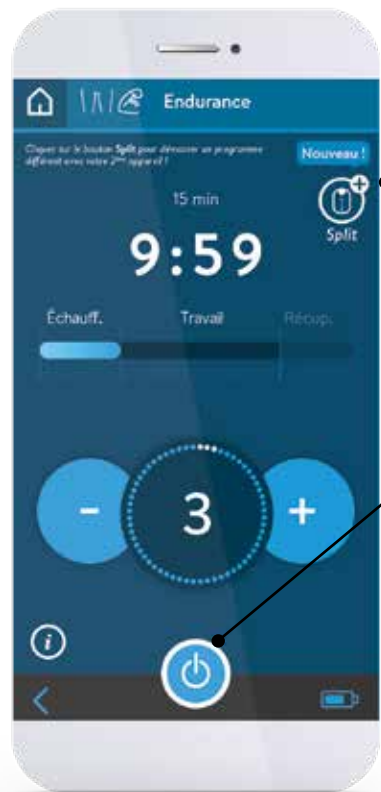
### DURING THE PROGRAM

> **Adjust the program's intensity turning Duo Sport wheel or using the App.**

To achieve greater effectiveness, follow our intensity advices.

### SOLO MODE / SPLIT FUNCTION

> **If you are using Duo Sport in SOLO mode and you wish to start a different program on the second device, click on SPLIT icon to select a second program to run on your second Duo Sport device.**

### PROGRAM TERMINATION

> **If you want to stop the program, press the «Stop» button.** You will be redirected to the main menu.

### END OF THE PROGRAM

> **A message will indicate the end of the program.** Turn Duo Sport off before removing the electrodes. Reposition the electrodes on the proper side of the plastic film and store them inside the sealed bag.

## SAFETY INSTRUCTIONS & TECHNICAL CHARACTERISTICS

### CONTRAINDICATIONS

The use of Duo Sport is forbidden for people fitted with a pacemaker or epilepsy-prone. Duo Sport use is also not recommended for pregnant women (⚠️), and people with a hernia or an incisional hernia. People suffering from heart weakness must avoid areas around the heart. Don't use Duo Sport if you have heart irregularities or heart diseases.

### ⚠️ IMPORTANT !

Do not use Duo Sport in water or in humid environments (sauna, balneotherapy and hammam, etc.). It is forbidden to place the electrodes on the sinocarotid region (in the front of the neck), on the transcranial region (through the head) or transthoraxal region (through the ribcage). To get it right, please follow the electrode placement's advices, they have been prepared by health professionals and correspond to areas of the body where Duo Sport will be more effective.

- Weight: 25 grams
- Battery : 3,7V 400mAh
- Voltage : 0-70V (500Ω load)
- Frequency : 1Hz-1200Hz
- Output Current : Maximum 300mA (pic), average of 20 mA
- Output voltage : Maximum 120V (500Ω load), 65V with load
- Power supply : USB
- USB (on charge) : 5V 500mA, 2.5W
- Impulse form : symmetrical biphasic square pulse
- Maximum resistance : 1500 Ohm
- Automatic shut-off after each program : yes

Information on waste management

At the end of the product life, make sure to discard it properly, following the relevant national rules. We thank you, in advance, for dropping the product in a place provided for this purpose in order to ensure a safe disposal in the respect of the environment.

### RECOMMENDATIONS (storage and transport)

Temperature : 10°C / 45°C - Humidity level : 25% / 80%

Duo Sport must be used only in accordance with the proposed programs. Bluetens Limited cannot be held responsible for any misuse or manipulation of Duo Sport. Handle electrodes only when Duo Sport is turned off. Duo Sport is a medical device, do not leave it near children. Your Duo Sport is a class 2A medical device certified by TÜV. It comes with a 2 year warranty. Hereby, Bluetens limited, declares that Duo Sport is in compliance with the essential requirements and other relevant provisions of 1999/5/EC directive and 3/42/EEC directive 9. A copy of the Declaration of conformity is available on request to Bluetens Ltd, 2003 20/F Tower 5, China Hong Kong City, 33 Canton Road, TST, Hong-Kong ou Bluetens France, 5 passage Saint Bernard, 75011 Paris, France and on [www.bluetens.com/doc](http://www.bluetens.com/doc)

### SAFETY PRECAUTIONS

Do not use Duo Sport within a meter from shortwave medical devices or microwave equipment, which could modify the current generated by the device. Use only the cable provided by Duo Sport. During sessions do not disconnect the electrode cable because power is on. The device must be turned off first. Never take off the electrodes when Duo Sport is powered on. Never connect the electrode cable to an external power source. That could cause electric shock. Do not use any damaged or defective electrode cable. Always use the USB cable provided by Duo Sport to charge the device. Never use Duo Sport or the USB cable if damaged. Keep Duo Sport and its accessories out of reach of children. Prevent any foreign matter (soil, water, metal, etc.) from entering the device. Do not use Duo Sport while driving. Never recharge Duo Sport during its use; a sliding part (the slider) prevents the simultaneity of recharge and use as an extra protection. In case of doubt, do not hesitate to consult a doctor.

- Ignition safety : 100% (⚠️)
- Classified as type BF (⚠️)
- CE regulations concerning Medical Devices (Directive 93/42/EEC) (CE 0123)
- ROHS : devices built without materials that can have a harmful impact on the environment. **RoHS**
- Protection against electric shock class 2 (⚠️)
- Device in contact with the patient out heart area (⚠️)
- Read the label (⚠️)

## CUSTOMER SERVICE

We manage all after-sale questions live.

If you need assistance after the purchase of your Duo Sport from one of our distributors, do not return the product to the physical store, please contact us directly at [contact@bluetens.com](mailto:contact@bluetens.com)

Our agents will get back to you within 24 hours with a quick solution.

## GUARANTEE

Duo Sport is guaranteed for two years, this warranty does not cover damages caused by improper maintenance or usage not compliant with the related documentation, accidents (shock, fire, liquid), modifications, repairs, misuse or normal wear and tear. This guarantee covers material or manufacturing defects. This guarantee covers products purchased directly from Bluetens or from a Bluetens authorized dealer only.

### GUARANTEE PERIOD

Your product is guaranteed two years from the date of purchase, the battery has a six month guarantee.

### NOT COVERED BY THE GUARANTEE

Electrodes, clip, satchel, as well as cables and accessories included in the box do not benefit from this guarantee.

### BLUETENS AT YOUR SERVICE

In the event of a defect covered by the guarantee, Bluetens will send you a return slip and will be responsible for the repairs or exchange of your device as well as of its return to you.

### HOW TO GET CUSTOMER SUPPORT?

You can reach our Customer Service at any time at the following address:  
[contact@bluetens.com](mailto:contact@bluetens.com), a member of our team will be dedicated to find a solution to your request.

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# User instructions for TENS

### CAUTION

Do not use while sleeping.  
Do not use while driving.  
Do not continuously keep electrodes on body.  
Do not use metal to touch the electrode.  
Use this product with high frequency surgery devices might cause those device to become unstable.

### MAIN DEVICE MAINTENANCE

- Clear regularly.
- Take off cables from the device and electrode after power off.
- Use soft cloth to clean if the device is dirty.
- Use wet cloth with neutral detergent or spirit to clean if necessary.
- Do not use spirit on the LED panel.
- Do not use other detergent.
- Keep water or liquids away from the device.

### OPERATION INSTRUCTIONS

The device needs to be used together with the water proof carrying case in order to prevent any damage caused by water. Please put the device into the carrying bag during use. Then use the App to realize all settings and controlling of the device.

## ELECTRODE PAD

Specification :

ITEM	STANDARD
Size	49.5±0.5 mm
Thickness	1.6 + 0.4 mm
Appearance	No stain, no gel spillover at the edges
Backing material	White non-woven cloth
Gel adhesion	Under normal humidity (60%), skin application is above 100 times after 72 hour unprotected storing
	Under high humidity (80%), skin application is above 30 times after 72-hour unprotected storing
	Under low humidity (20%), skin application is above 10 times after 72-hour unprotected storing
	Dried completely at 100°C for 1.5 hours, skin application is above 5 times
	For un-pretreated or sweaty skin, skin application is above 10 times.
	Comfortable cutaneous sensation
	No edge lift after 2-hour skin application under any condition
	Skin adhesion reaction time: No need to press or press for less than 3 seconds

ITEM	STANDARD
Conductivity	After 100 skin application, peel the gel off from the carbon film in the water and dry it. The impedance is <200ohms between any two points
	A.C. impedance for gel-to-metal by LCR meter: <100ohms at 100Hz; <80ohms at 1,000Hz <70ohms at 10,000Hz; <65ohms at 100,000Hz
Lead wire	Adaptable to male pin from 1.95~2.10m
	Strength between leadwire and connector is min. 4kg pulling force
	Strength between leadwire and electrode pad is min.3kg pulling force
Adhesive force	Gel is strongly attached to the carbon film, hard to peel it off
Protective liner	110X110mm, 0.075mm transparent release liner, no gel residue when peel the electrode off
Printing	Lot number & expiry date: firm and legible
Safety	Biocompatibility of the gel G607 used for this product conforms to ISO10993-5-10
Package	120±5mmX185±5mm, transparent eco-friendly PE re-sealable bag
Applicable humidity	Humidity: 15%~70% Temperature: 0°C~30°C
Shelf life	2-year sealed storage under 30°C, gel color and adhesion remain unchanged

## ELECTROMAGNETIC COMPATIBILITY DESCRIPTIONS

The EMC declaration according to the requirement of EN 60601-1-2

### CAUTIONS

User must regard EMC, please install and put in service BLT15DS according to the EMC information provided in the accompanying documents.

Portable and mobile RF communications equipment can affect medical electrical equipment.

The performance of the EQUIPMENT and SYSTEM that was determined to be essential performance.

Table 201-Guidance and manufacturer's declaration-electromagnetic emissions for BLT15DS, as following table.

Table 202-Guidance and manufacturer's declaration-electromagnetic immunity for BLT15DS, as following table.

Table 204 -Guidance and manufacturer's declaration-electromagnetic immunity for BLT15DS, as following table.

Table 206 Recommended separation distances between portable and mobile RF communications equipment and BLT15DS, as following table.

Table 201

GUIDANCE AND MANUFACTURER'S DECLARATION – ELECTROMAGNETIC EMISSIONS		
The BLT15DS is intended for use in the electromagnetic environment specified below. The customer or the user of the BLT15DS should assure that it is used in such an environment.		
Immunity test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The BLT15DS uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The BLT15DS is suitable for use in all establishments, including domestic establishments.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Not applicable	

Table 204

GUIDANCE AND MANUFACTURER'S DECLARATION – ELECTROMAGNETIC IMMUNITY			
The BLT15DS is intended for use in the electromagnetic environment specified below. The customer or the user of the BLT15DS should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	Not applicable	Portable and mobile RF communications equipment should be used no closer to any part of the BLT15DS, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance 80 MHz to 800 MHz 800 MHz to 2,5 GHz  Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, a should be less than the compliance level in each frequency range's  Interference may occur in the vicinity of equipment marked with the following symbol:  
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the BLT15DS is used exceeds the applicable RF compliance level above, the BLT15DS should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the BLT15DS. b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Table 206

RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT AND THE BLT15DS			
The BLT15DS is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the BLT15DS can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the BLT15DS as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter  W	Separation distance according to frequency of transmitter		
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2,5 GHz
	$d = 1,2\sqrt{P}$	$d = 1,2\sqrt{P}$	$d = 2,3\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance <i>d</i> in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where <i>P</i> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

[contact@bluetens.com](mailto:contact@bluetens.com)  
[www.bluetens.com](http://www.bluetens.com)